

## FFM-N8040GA-1K

FFM Acid-Resistant Nanofiltration Membranes selective separation for acid whey demineralization, citrus juice decolorization, antibiotic polishing, and nickel/cobalt recovery. Features with pH 0-12 stability, 1,200+ acid cycles, >90% divalent ion rejection retention, and minimal MWCO drift-ensuring compliance and yield in aggressive acidic processes.

### Performance Parameters


Membrane Material:	Special Composite Materials
Structure:	Sanitary Spiral-Wound Mesh Design
Effective Retention Accuracy:	1000 Dalton
Permeate Tube:	CPVC
Standard Treatment Solution:	10% HCL/10% H <sub>2</sub> SO <sub>4</sub> /15% H <sub>3</sub> PO <sub>4</sub> /2% HNO <sub>3</sub>

Item	31 mil	46 mil
Maximum Salt rejection (MgSO <sub>4</sub> )	97.0%	97.0%
Minimum Salt rejection (MgSO <sub>4</sub> )	95.0%	95.0%
Standard Permeate Flow	6500 gpd (24.6 m <sup>3</sup> /d)	5000 gpd (18.9 m <sup>3</sup> /d)

Test conditions:	2000 mg/L MgSO <sub>4</sub> Solution
	86 F(30 °C) Operating Temperature
	145 psi (1.0 Mpa) Applied Pressure
	15% Permeate Recovery

\*Single-element water outlet tolerance: ±20% (standard conditions).

### Product Specifications

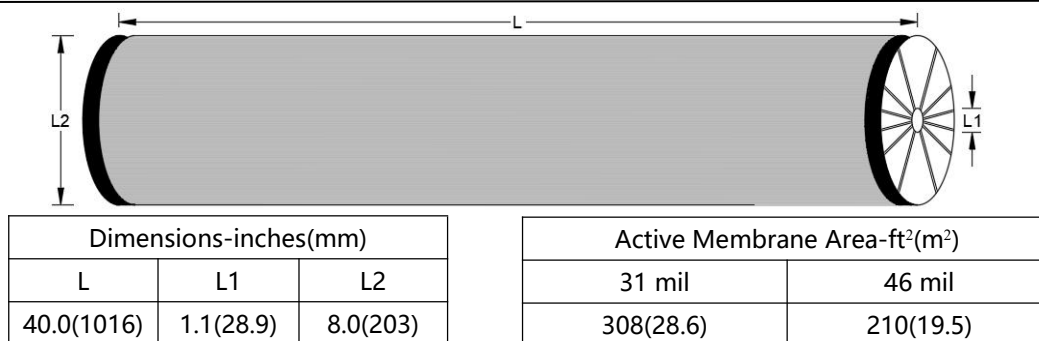

### Operation Parameters

Maximum Operating Temperature:	122 °F (50 °C)
Recommended Operating Pressure:	220-510psi (1.5-3.5 Mpa)
pH Range,Continuous Operation:	0.0-12.0(25°C)
pH Range, Ceb:	0.0-12.0(25°C)
Design Pressure Drop per Membrane Element:	10 psi (0.07 Mpa)
Design Pressure Drop per Membrane Housing (5 Elements in Series):	60 psi (0.42 Mpa)

\*For specialized applications or customized solutions, please contact FFM Inc. directly.